

Trypsin-EDTA (0.05 %) in DPBS (1x) LTDe02

General Information

Trypsin-EDTA solutions are used to detach adherent cells from culture surfaces. They are composed of natural porcine pancreas-derived trypsin and EDTA. The concentration of trypsin necessary to dislodge cells from their substrate is dependent primarily on the cell type and the age of the culture. Various formulations should be tested to determine the best product for a specific application.

Appearance	Clear frozen liquid
Storage and shelf life	Store at $\leq -15^{\circ}\text{C}$. Avoid repeated freeze-thaw cycles. Preparation of aliquots recommended. Once opened, store at $+4^{\circ}\text{C}$ and use within 2-4 weeks.
Shipping conditions	Frozen (Dry ice)
Thawing	$+37^{\circ}\text{C}$ water bath or overnight at $+2^{\circ}\text{C}$ to $+8^{\circ}\text{C}$. Swirl gently to homogenize.

Formulation


Components	mg/l
EDTA 4Na	220
KCl	200
KH_2PO_4	200
NaCl	8000
Na_2HPO_4	1150
Trypsin	500

Instructions for Use

Detachment of adherent cells using Trypsin-EDTA

Trypsin-EDTA (0.05%) in DPBS (1x) solution is supplied as a sterile, ready-to-use, frozen liquid. This entire procedure should be done in a laminar flow hood using proper aseptic technique.

1. Carefully aspirate all of the media from the cell culture flask.
2. Rinse cells with Ca^{2+} and Mg^{2+} -free salt solution (see related products), aspirate, and discard.
3. Pre-warm the trypsin solution in a $+37^{\circ}\text{C}$ water bath. Add enough trypsin solution to completely cover the cells.
4. Incubate the flask at $+37^{\circ}\text{C}$, or for more sensitive cultures, at room temperature or $+2^{\circ}\text{C}$ to $+8^{\circ}\text{C}$.
5. When the trypsinization process is complete, cells will appear rounded upon microscopic examination and the solution in the flask will appear cloudy. Check the flask often to avoid overexposure. Trypsin can cause cellular damage and time of exposure should be kept to a minimum.
6. The time required to detach cells from the culture surface is dependent on the cell type, the age of the culture, population density, serum concentration in the growth medium and time since last subculture.

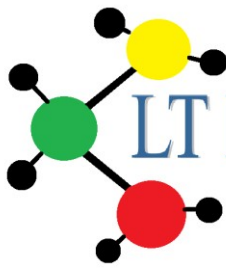

UAB "LT Biotech"

Į.K. 302303586, PVM kodas LT 100004741118,
Rugių 21-24, LT-08419, Vilnius Reg. Nr. 127918,
V.Į. Registrų centras Vilniaus filialas
Tel/fax +370 5 216 02 27

7. Neutralizetrypsineitherwithserumcontainingmediumortrypsininhibitor.Gentlycentrifugethecellsuspension and discard the trypsin-containing supernatant.
8. Resuspendthecellpelletwithfreshmediumandcountorcultureasdesired.

Precautions and Disclaimer

This product is for research use only.